HAEMATOKRIT 210

ⒹⒺ	Bedienungsanleitung	10
Ⓔ⒩	Operating Instructions.....	25
ⒻⓇ	Mode d'emploi	40
ⒾⓉ	Istruzioni per l'uso	55

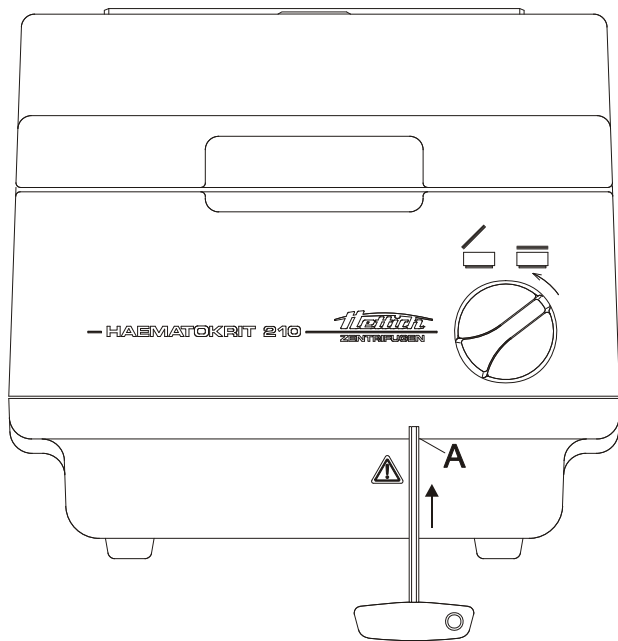
Fig. 1

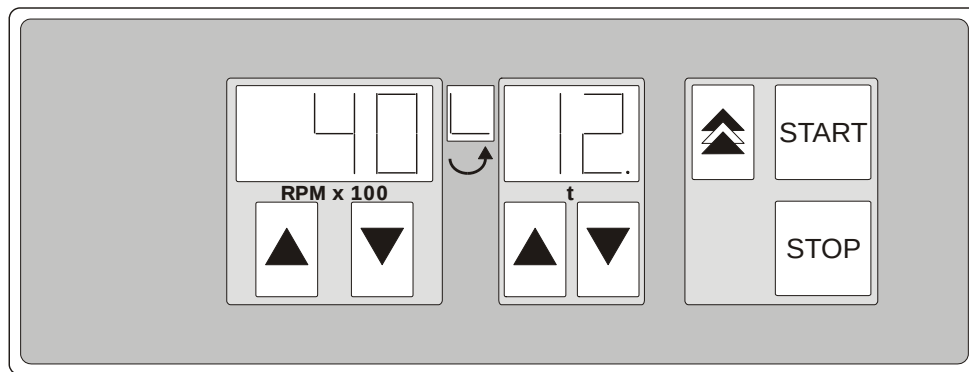
Fig. 2

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EC Declaration of conformity

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des Herstellers / of the manufacturer / du fabricant / del costruttore
Andreas Hettich GmbH & Co. KG • Föhrenstraße 12 • D-78532 Tuttlingen • Germany

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We hereby declare under our sole responsibility that the designated device and its accessories, which are listed in the technical documentation for this device and whose conformity has been assessed together with the device, conform to the Directive 98/79/EC on in vitro diagnostic medical devices.

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Device for centrifugal separation of sample material of human origin for preparation of in-vitro diagnostic analyses.

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The conformity evaluation process was performed in accordance with appendix III of Directive 98/79/EC.

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Applied standards and directives:

According to the list of applied standards and valid directives which is part of the product documentation.

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Tuttlingen, 2012-05-03

A handwritten signature in black ink, appearing to read 'H. Eberle', written in a cursive style.

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- **USA:** QSR, 21CFR 820 "CFR Title 21 - Food and Drugs: TITLE 21- FOOD AND DRUGS, CHAPTER I - FOOD AND DRUG ADMINISTRATION DEPARTMENT OF HEALTH AND HUMAN SERVICES, SUBCHAPTER H - MEDICAL DEVICES, Part 820 QUALITY SYSTEM REGULATIONS"
- **Kanada:** CMDR, SOR/98-282 "Medical Devices Regulations"

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Standards and regulations which apply to this device

The device is a high-end technical product. It is subject to extensive testing and certification procedures according to the following standards and regulations in their respectively valid version:

Electrical and mechanical safety for design and final testing:

Standard series: IEC 61010 (conform to standards of DIN EN 61010)

- IEC 61010-1 "Safety requirements for electrical equipment for measurement, control, and laboratory use - Part 1: General requirements" (Pollution Degree 2, Installation Category II)
- IEC 61010-2-010 "Safety requirements for electrical equipment for measurement, control and laboratory use - Part 2-010: Particular requirements for laboratory equipment for the heating of materials" (applied to heated centrifuges only)
- IEC 61010-2-020 "Safety requirements for electrical equipment for measurement, control, and laboratory use - Part 2-020: Particular requirements for laboratory centrifuges"
- IEC 61010-2-101 "Safety requirements for electrical equipment for measurement, control and laboratory use - Part 2-101: Particular requirements for in vitro diagnostic (IVD) medical equipment"

Electromagnetic Compatibility:

- EN 61326-1 "Electrical equipment for measurement, control and laboratory use - EMC requirements - Part 1: General requirements"

The EMC standard refers to the following generic standards:

Emission:

- | | |
|------------------------------------|------------------|
| • Conducted emission | EN 55011 Class B |
| • Radiated emission | EN 55011 Class B |
| • Harmonic current emissions | EN 61000-3-2 |
| • Voltage fluctuations and flicker | EN 61000-3-3 |

Immunity to...

- | | |
|--|---------------|
| • Electrostatic discharge | EN 61000-4-2 |
| • Electromagnetic fields | EN 61000-4-3 |
| • Electrical fast transients (Burst) | EN 61000-4-4 |
| • Surge | EN 61000-4-5 |
| • Conducted RF disturbances | EN 61000-4-6 |
| • Voltage dips and short interruptions | EN 61000-4-11 |

European directives applied for conformity assessment procedures:

In vitro diagnostic device directive 98/79/EG

EC conformity assessment procedure according to annex III "EC DECLARATION OF CONFORMITY" – self-declaration by the manufacturer

Further partly applicable European directives:

- Machinery Directive 2006/42/EC
- EMC directive 2004/108/EC
- Low voltage directive 2006/95/EC

Applied medical device regulations outside Europe:

- **USA:** QSR, 21CFR 820 "CFR Title 21 - Food and Drugs: TITLE 21- FOOD AND DRUGS, CHAPTER I - FOOD AND DRUG ADMINISTRATION DEPARTMENT OF HEALTH AND HUMAN SERVICES, SUBCHAPTER H - MEDICAL DEVICES, Part 820 QUALITY SYSTEM REGULATIONS"
- **Canada:** CMDR, SOR/98-282 "Medical Devices Regulations"

Certified quality management system according to

- ISO 9001 "Quality management systems – Requirements"
- ISO13485 "Medical devices - Quality management systems - Requirements for regulatory purposes"

Environmental management system according to

- ISO 14001 "Environmental management systems - Requirements with guidance for use"

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Contents

1	Use according to specification	26
2	Remaining risks.....	26
3	Technical specifications	26
4	Notes on safety	27
5	Symbol meanings	29
6	Delivery checklist	29
7	Initial operation.....	30
8	Opening and closing the lid.....	30
8.1	Opening the lid.....	30
8.2	Closing the lid	30
9	Control and display elements.....	31
9.1	Symbols on the control panel	31
9.2	Keys and setting options.....	31
10	Setting the brake step	31
11	Centrifugation.....	32
11.1	Centrifugation with preselected time.....	32
11.2	Continuous operation.....	32
11.3	Short-time centrifugation.....	32
12	Relative centrifugal force (RCF).....	33
13	Centrifugation of materials or mixtures of materials with a density higher than 1.2 kg/dm ³	33
14	Haematocrit measurement.....	33
14.1	Haematocrit capillaries - Overview	33
14.2	Filling of the haematocrit capillaries.....	34
14.3	Sealing and centrifugation of the haematocrit capillaries	34
14.3.1	Standard capillary tubes, Order No.: 2074, and mylar-reinforced capillary tubes, Order No.: 1072	34
14.3.2	Self-sealing capillary tubes (Order No. 1071).....	34
14.4	Loading the rotor	34
14.5	Determination of the centrifugation time	35
14.6	Reading the haematocrit value	35
15	Emergency release	35
16	Maintenance and servicing	35
16.1	Centrifuge (housing, lid and centrifuging chamber)	36
16.1.1	Surface cleaning and care	36
16.1.2	Surface disinfection	36
16.1.3	Removal of radioactive contaminants.....	36
16.2	Cleaning and disinfecting the haematocrit rotor	36
16.3	Procedure in case of broken haematocrit capillaries	36
16.4	Rotor removal / installation	37
17	Faults	38
18	Change mains input fuses.....	39
19	Returning Devices.....	39
20	Disposal	39
21	Anhang / Appendix.....	70
21.1	Rotoren und Zubehör / Rotors and accessories	70

1 Use according to specification

This device is a medical product (laboratory centrifuge) within the context of the IVD Directive 98/79/EC.

The centrifuge is used for centrifugal separation of human blood samples in haematocrit capillaries in accordance with EN ISO 12772 or in special manufacturer-specific capillaries for quantitative Buffy Coat Analysis.

The centrifuge is only meant for this purpose.

Another use or one which goes beyond this, is considered to be non-intended. The company Andreas Hettich GmbH & Co. KG is not liable for damage resulting from this.

Observing all information in the operating instructions and complying with the measures described therein is also a part of the intended use.

2 Remaining risks

The device is built according to the state-of-the-art and the recognized safety regulations. If used and handled improperly, there could be life-threatening danger to the user or third parties, or the device could be impaired or there could be other property damage. The device is only to be used for its intended purpose and only when it is in safe working condition.

Malfunctions which could affect safety must be corrected immediately.

3 Technical specifications

Manufacturer	Andreas Hettich GmbH & Co. KG D-78532 Tuttlingen	
Model	HAEMATOKRIT 210	
Type	2104	2104-01
Mains voltage (± 10%)	208 - 240 V 1~	100 - 127 V 1~
Mains frequency	50 - 60 Hz	50 - 60 Hz
Connected load	250 VA	220 VA
Current consumption	1.35 A	2.2 A
Max. capacity	24 x 2.2 ml	
Allowed density	1.2 kg/dm ³	
Speed (RPM)	13000	
Force (RCF)	16060	
Kinetic energy	850 Nm	
Obligatory inspection (BGR 500)	no	
Ambient conditions (EN / IEC 61010-1) – Set-up site – Altitude – Ambient temperature – Humidity – Excess-voltage category (IEC 60364-4-443) – Pollution degree	Indoors only Up to 2000 m above sea level 2°C to 40°C Maximum relative humidity 80% for temperatures up to 31°C, linearly decreasing to 50% relative humidity at 40°C. II 2	
Device protection class	I	
Not suitable for use in explosion-endangered areas.		
EMC – Emitted interference, Interference immunity	EN / IEC 61326-1, Class B	FCC Class B
Noise level (dependent on rotor)	≤ 58 dB(A)	
Dimensions – Width – Depth – Height	275 mm 330 mm 247 mm	
Weight	approx. 10 kg	

4 Notes on safety

No claim of warranty will be considered by the manufacturer unless ALL instructions in this manual have been followed.

- The centrifuge should be installed on a good, stable base.
- Before using the centrifuge absolutely check the rotor for firm placement.
- When the centrifuge is running, according to EN / IEC 61010-2-020, no persons, dangerous substances or objects may be within the safety margin of 300 mm around the centrifuge.
- Rotors, suspensions and accessories that possess traces of corrosion or mechanical damage or if their term of use has expired may not be used any longer.
- The centrifuge may no longer be put into operation when the centrifuging chamber has safety-related damages.
- With swing-out rotors the trunnions must be regularly lubricated (Hettich Lubricating Grease No. 4051) in order to ensure consistent swinging out of the hangers.
- For centrifuges without temperature control, when the room temperature is increased and/or if the device is frequently used, the centrifuging chamber could be heated up. Therefore, it can't be ruled out that the sample material might be changed due to the temperature.

- Before the initial operation of your centrifuge you should read and pay attention to the operating instructions. Only personnel that has read and understood the operating instructions are allowed to operate the device.
- Along with the operating instructions and the legal regulations on accident prevention, you should also follow the recognised professional regulations for working in a safe and professional manner. These operating instructions should be read in conjunction with any other instructions concerning accident prevention and environmental protection based on the national regulations of the country where the device is to be used.
- This centrifuge is a state-of-the-art piece of equipment which is extremely safe to operate. However, it can lead to danger for users or others if used by untrained staff, in an inappropriate way or for a purpose other than that it was designed for.
- The centrifuge must not be moved or knocked during operation.
- In case of fault or emergency release, never touch the rotor before it has stopped turning.
- To avoid damage due to condensate, when changing from a cold to a warm room the centrifuge must either heat up for at least 3 hours in the warm room before being connected to the mains, or run hot for 30 minutes in the cold room.
- Only the rotors and accessories approved by the manufacturer for this device may be used (see chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories"). Before centrifuge vessels are used which are not listed in the chapter "Appendix, Rotors and accessories", the user must make sure they can be used by asking the manufacturer.
- The centrifuge rotor may only be loaded in accordance with the chapter "Loading the rotor".
- When centrifuging with maxim revolutions per minute the density of the materials or the material mixtures may not exceed 1.2 kg/dm³.
- The centrifuge may only be operated when the balance is within the bounds of acceptability.
- The centrifuge may not be operated in explosion-endangered areas.
- The centrifuge must not be used with:
 - inflammable or explosive materials
 - materials that react with one another producing a lot of energy.

- If users have to centrifuge hazardous materials or compounds contaminated with toxic, radioactive or pathogenic micro-organisms, they must take appropriate measures.
For hazardous substances centrifuge containers with special screw caps must strictly be used. In addition to the screw cap centrifuge containers, for materials in hazard category 3 and 4 a biosafety system must be used (see the World Health Organisation's "Laboratory Biosafety Manual").
In a biosafety system, droplets and aerosols are prevented from escaping by a bioseal (packing ring).
If the hanger of a biosafety system is used without the lid, the packing ring must be removed from the hanger in order to prevent the packing ring from being damaged during the centrifugation run. Damaged packing rings must not be used to seal the biosafety system.
Without the use of a biosafety system the centrifuge is not microbiologically sealed in the sense of the EN / IEC 610101-2-020 standard.
When closing a biosafety system, please observe the instructions provided in the chapter, "Closing Biosafety Systems".
For further details of available biosafety systems see chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories". If in doubt, you should obtain relevant information from the manufacturer.
- The centrifuge must not be operated with highly corrosive substances which could impair the mechanical integrity of rotors, hangers and accessories.
- Repairs must only be carried out by personnel authorised to do so by the manufacturer.
- Only original spare parts and original accessories licensed by the Andreas Hettich GmbH & Co. KG company are allowed to be utilised.
- The following safety regulations apply:
EN / IEC 61010-1 and EN / IEC 61010-2-020 as well as their national deviations.
- The safe operation and reliability of the centrifuge can only be guaranteed if:
 - the centrifuge is operated in accordance with the operating instructions,
 - the electrical installation on the site where the centrifuge is installed conforms to the demands of EN / IEC stipulations,
 - the tests for device safety required in the respective countries, e.g. in Germany in acc. with BGV A1 and BGR 500, are carried out by an expert.

5 Symbol meanings

Symbol on the device:

Attention, general hazard area.

Before using the device, make sure you read the operating instructions and observe the safety information!

Symbol in this document:

Attention, general hazard area.

This symbol refers to safety relevant warnings and indicates possibly dangerous situations.

The non-adherence to these warnings can lead to material damage and injury to personal.

Symbol on the device and in this document:

Beware of biohazard.

Symbol in this document:

This symbol refers to important circumstances.

Symbol on the device and in this document:

Symbol for the separate collection of electric and electronic devices according to the guideline 2002/96/EG (WEEE). The device belongs to Group 8 (medical devices).

Applies in the countries of the European Union, as well as in Norway and Switzerland.

6 Delivery checklist

The following items and accessories are delivered with the centrifuge:

- 1 Connecting cable
- 2 Fuses
- 1 Notes on moving the equipment safely
- 1 Operating instructions
- 1 Hexagonal pin-type spanner

The rotor(s) and associated accessories are included in the delivery in the quantity.

7 Initial operation

- Remove the transportation safety device from the bottom of the housing, see sheet "Transportation safety device"
- **Position the centrifuge in a stable and level manner in a suitable place. During set-up, the required safety margin of 300 mm around the centrifuge is to be kept according to EN / IEC 61010-2-020.**

When the centrifuge is running, according to EN / IEC 61010-2-020, no persons, dangerous substances or objects may be within the safety margin of 300 mm around the centrifuge. Substances can be discharged through the ventilation opening on the rear side of the centrifuge. The device is to be placed in such a way that the air stream is not directed towards people.

- Do not place any object in front of the ventiduct.
Keep a ventilation area of 300 mm around the ventiduct.
- Check whether the mains voltage tallies with the statement on the type plate.
- Connect the centrifuge with the connection cable to a standard mains socket. For connection ratings refer to Chapter "Technical specifications".
- Turn on the mains switch. Switch position "I".
The last used centrifuge data will be displayed.
- Open the lid.
- Remove the transportation safety device in the centrifuge, see sheet "Transportation safety device".

8 Opening and closing the lid

8.1 Opening the lid

The lid can only be opened when the centrifuge is switched on and the rotor is at rest. If it cannot be opened under these circumstances, see the section on "Emergency release".

- Rotate the turning handle on the front panel to the left. The symbol "L" (lid open) illuminates in the rotation indicator ↶.
- Open the lid.

8.2 Closing the lid

Do not bang the lid shut.

- Place the lid and lightly press down the front edge of the lid.
- Turn the handle on the front panel to the right. The symbol "—" (lid closed) illuminates in the rotation indicator ↷.

9 Control and display elements

See figure on page 2.

Fig. 2: Display and control panel

9.1 Symbols on the control panel

Rotation indicator. The rotation indicator lights up and rotates anticlockwise while the rotor is turning.

When the rotor is stationary, the status of the lid is displayed by symbols in the rotation indicator:

Symbol : Lid open

Symbol : Lid closed

Operator errors and occurring faults are indicated on the display (see Chapter "Faults").

9.2 Keys and setting options

RPM x 100

- Speed
A numeric value of 500 RPM up to the maximum rotor speed can be set. For maximum rotor speed, see chapter "Anhang/Appendix, Rotor und Zubehör/Rotor and accessories". Preset in steps of 100 (RPM = displayed value x 100).
If the key or is kept pressed, the value changes with increasing speed.
- Display the brake step.

t

- Running time
- Preset from 1 - 99 minutes, in 1 minute steps
- Continuous operation "--"
• Braking steps 0 or 1. Step 1 = short run-down time, Step 0 = long run-down time.
If the key or is kept pressed, the value changes with increasing speed.

- Start centrifugation run.


- End centrifugation run.
The rotor runs down with the preselected brake step.
- Save the brake step.

- Short-time centrifugation.
The centrifugation run occurs while the key is kept pressed.
- Display the brake step.

10 Setting the brake step

- Switch off the mains switch.
- Keep the key beneath the speed indicator and the key pressed simultaneously.
- Switch on the mains switch and release the keys again.
The speed indicator shows the machine version and the time indicator shows the set brake step: e.g.:

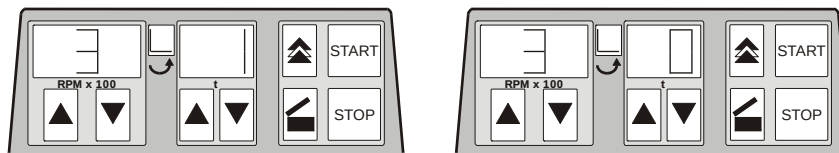If the machine version and brake step are not displayed, press the key under the speed indicator until they are displayed.

The machine version is set by the manufacturer and cannot be changed.

- Set the desired brake step with the keys beneath the time indicator.
Step 1 = short run-down time, Step 0 = long run-down time.
For run-down times, see chapter "Anhang/Appendix, Rotor und Zubehör/Rotor and accessories".
- Press the key to save the setting.

11 Centrifugation

When the centrifuge is running, according to EN / IEC 61010-2-020, no persons, dangerous substances or objects may be within the safety margin of 300 mm around the centrifuge.

If the permissible weight difference is exceeded within the rotor loading, the drive switches off during the run-up time, and error -3- is displayed (see chapter "Faults").

The centrifugation run can be interrupted at any time by pressing the key **STOP**.

The time and speed can be changed during the centrifugation run, with the keys **▲** **▼**.

If the key **▲** or **▼** is kept pressed, the value changes with increasing speed.

After a centrifugation run, the display flashes until the cover is opened or a key is pressed.

If the symbol "—" (lid closed) and "L" (lid open) flashes alternately in the rotation indicator , operation of the centrifuge can only be continued after opening the lid.

- Switch on the mains switch (switch position "I").
- Load the rotor and close the centrifuge cover.

11.1 Centrifugation with preselected time

- Set the desired speed with the keys **▲** **▼** beneath the speed indicator.
- Set the desired time with the keys **▲** **▼** beneath the time indicator.
- Press the key **START**. The rotation indicator  appears while the rotor is turning.

The time is displayed in minutes. The last minute is counted down in seconds. When the time is displayed in minutes, a point flashes next to the number.

- After expiry of the time or if the centrifugation run is interrupted by pressing the key **STOP**, the rotor runs down with the set brake step.

During the centrifugation run the rotor speed and the remaining time are displayed.

11.2 Continuous operation

- Set the desired speed with the keys **▲** **▼** beneath the speed indicator.
- Set the time to zero with the key **▼** beneath the time indicator. "--" is displayed.
- Press the key **START**. The rotation indicator  appears while the rotor is turning. The time count starts from 0.

The first minute is counted up in seconds, and then the time is displayed in minutes. When the time is displayed in minutes, a point flashes next to the number.

- Press the key **STOP** to end the centrifugation run. The rotor runs down with the set brake step.

During the centrifugation run the rotor speed and the expired time are displayed.

11.3 Short-time centrifugation

- Set the desired speed with the keys **▲** **▼** beneath the speed indicator.
- Keep the key  pressed. The rotation indicator  appears while the rotor is turning. The time count starts from 0.

The first minute is counted up in seconds, and then the time is displayed in minutes. When the time is displayed in minutes, a point flashes next to the number.

- Release the key  again to end the centrifugation run. The rotor runs down with the set brake step.

During the centrifugation run the rotor speed and the expired time are displayed.

12 Relative centrifugal force (RCF)

The relative centrifugal force (RCF) is given as a multiple of the acceleration of gravity (g). It is a unit-free value and serves to compare the separation and sedimentation performance.

These values are calculated using the formula below:

$$\text{RCF} = \left(\frac{\text{RPM}}{1000} \right)^2 \times r \times 1,118 \quad \Rightarrow \quad \text{RPM} = \sqrt{\frac{\text{RCF}}{r \times 1,118}} \times 1000$$

RCF = relative centrifugal force

RPM = rotational speed (revolutions per minute)

r = centrifugal radius in mm = distance from the centre of the turning axis to the bottom of the centrifuge.
For more on the centrifugal radius see the chapter "Anhang/Appendix, Rotoren und Zubehör/Rotors and accessories".

The relative centrifugal force (RCF) stands in relation to the revolutions per minute and the centrifugal radius.

13 Centrifugation of materials or mixtures of materials with a density higher than 1.2 kg/dm³

When centrifuging with maximum revolutions per minute the density of the materials or the material mixtures may not exceed 1.2 kg/dm³.

The speed must be reduced for materials or mixtures of materials with a higher density.

The permissible speed can be calculated using the following formula:

$$\text{Reduced speed (n}_{\text{red}}) = \sqrt{\frac{1.2}{\text{Greater density [kg/dm}^3]}} \times \text{maximum speed [RPM]}$$

e.g.: maximum speed RPM 4000, density 1.6 kg/dm³

$$n_{\text{red}} = \sqrt{\frac{1.2 \text{ kg/dm}^3}{1.6 \text{ kg/dm}^3}} \times 4000 \text{ RPM} = 3464 \text{ RPM}$$

If in doubt you should obtain clarification from the manufacturer.

14 Haematocrit measurement

14.1 Haematocrit capillaries - Overview

Standard capillary tubes, Order No.: 2074

The classical heparinized haematocrit capillary tube made of glass. Tried and tested for years.

Mylar-reinforced capillary tubes, Order No.: 1072

The heparinized glass capillary is coated with plastic. Should such a capillary break the glass splinters are retained in this coating. This means that the risk of injury for the user is reduced.

Self-sealing capillary tubes, Order No. 1071

These capillary tubes, which are also heparinized, no longer need to be sealed. They are already equipped with a seal. At first, the seal material is air-permeable so that the capillary can be filled. The seal will only become impermeable when it comes into contact with blood. Therefore, there will be no leaks during centrifugation.

14.2 Filling of the haematocrit capillaries

The capillaries have different volumes according to type.

The following are the exact data:

Capillary:	Standard	Mylar-reinforced	Self-sealing
Order No.:	2074	1072	1071
Total volume without seal (µl):	49.6	74.9	74.9
Internal diameter (mm):	0.92	1.12	1.12
Overall length (mm)	75	75	75

Therefore, more blood is needed for the capillaries No. 1072 and 1071 to be filled. The dimensions, however, are still in compliance with the requirements of DIN 12846 (standard for haematocrit capillaries) and the accuracy of the haematocrit values is not affected.

For the haematocrit measurement the capillaries are always filled up to $\approx 75\%$:

Capillary:	Standard	Mylar-reinforced	Self-sealing
Order No.:	2074	1072	1071
Volume when filled up to 75%(µl)	37.2	56.2	56.2

14.3 Sealing and centrifugation of the haematocrit capillaries

14.3.1 Standard capillary tubes, Order No.: 2074, and mylar-reinforced capillary tubes, Order No.: 1072

1. When filling the capillary ensure that the end of the capillary opposite the filler opening remains dry.
2. Seal at the dry end of the capillary:
Insert the capillary tube vertically into the seal up to the point where the end of the capillary touches the bottom of the seal plate.
Slightly tilt the capillary tube to the side and remove it from the seal.
3. Place the capillary tubes with the lute side facing outwards in the haematocrit rotor and apply the rotor lid and click into place.
4. Centrifuge.
5. Determine the haematocrit value.

14.3.2 Self-sealing capillary tubes (Order No. 1071)

1. Fill the capillary tube.
2. Let the column of blood settle in the direction of the seal.

The blood must touch the seal!

3. Tap the seal end of the capillary three times on the laboratory table!
This serves to intensify the contact with the seal material and to enhance the sealing process. Only in this manner the capillaries will be sealed off!
4. Place the capillary tubes with the lute side facing outwards in the haematocrit rotor and apply the rotor lid and click into place.
5. Centrifuge.
6. Determine the haematocrit value.

14.4 Loading the rotor

- Load opposite positions equally.
- Put lid onto rotor and click into place.
- Place the lid and lightly press down the front edge of the lid.
- Turn the handle on the front panel to the right.
If the symbol  is illuminated in the rotation display, the lid has been closed properly.

14.5 Determination of the centrifugation time

The following approximation formula is used for the determination of the centrifugation time:

$$\frac{100000}{\text{RCF}} = \text{time [min]}$$

Example:

Given an RCF of 16060 $\frac{100000}{16060} \approx 7 \text{ min}$ of centrifugation time.

For the calculation of the RCF see chapter "Relative Centrifugal Force (RCF)".

14.6 Reading the haematocrit value

Perform the set-up as follows:

1. hold the haematocrit rotor.
2. move the evaluation disk (lid of the haematocrit rotor) until the 100%-line of the evaluation disk is in line with the upper margin of the liquid column.
3. Read the haematocrit value at the upper margin of the erythrocyte column.

Repeat steps 1 to 3 until all capillaries have been analysed.

15 Emergency release

The lid cannot be opened in the event of a power failure. An emergency release has to be executed by hand.

For emergency release disconnect the centrifuge from the mains.
Open the lid only during rotor standstill.

See figure on page 2.

Switch off the mains switch (switch position "0").

Look through the window in the lid to be sure that the rotor has come to a standstill.

Position the centrifuge on the table sufficiently far forward so that the hole (Fig. 1, A) for the emergency release on the base of the casing can be reached.

- Insert the hexagonal pin driver into this hole, push it upwards and at the same time rotate the lid-lock handle to the left.
- Open the lid.

16 Maintenance and servicing

The device can be contaminated.

Pull the mains plug before cleaning.

Before any other cleaning or decontamination process other than that recommended by the manufacturer is applied, the user has to check with the manufacturer that the planned process does not damage the device.

- Centrifuges, rotors and accessories must not be cleaned in rinsing machines.
- They may only be cleaned by hand and disinfected with liquids.
- The water temperature must be between 20 – 25°C.
- Only detergents/disinfectants may be used which:
 - have a pH between 5 - 8
 - do not contain caustic alkalis, peroxides, chlorine compounds, acids and alkaline solutions
- In order to prevent appearances of corrosion through cleaning agents or disinfectants, the application guide from the manufacturer of the cleaning agent or disinfectant are absolutely to be heeded.

Haematocrit rotor, rotor lid, padded shells, and sealing rings are not autoclavable.

16.1 Centrifuge (housing, lid and centrifuging chamber)

16.1.1 Surface cleaning and care

- Clean the centrifuge housing and the centrifuging chamber regularly, using soap or a mild detergent and a damp cloth if required. For one thing, this services purposes of hygiene, and it also prevents corrosion through adhering impurities.
- Ingredients of suitable detergents: soap, anionic tensides, non-ionic tensides.
- After using detergents, remove the detergent residue by wiping with a damp cloth.
- The surfaces must be dried immediately after cleaning.
- In the event of condensation water formation, dry the centrifugal chamber by wiping out with an absorbent cloth.
- Lightly rub the rubber seal of the centrifuge chamber with talcum powder or a rubber care product after each cleaning.
- The centrifuging chamber is to be checked for damage once a year.

If damage is found which is relevant to safety, the centrifuge may no longer be put into operation. In this case, notify Customer Service.

16.1.2 Surface disinfection

- If infectious materials penetrates into the centrifugal chamber this is to be disinfected immediately.
- Ingredients of suitable disinfectants: ethanol, n-propanol, isopropyl alcohol, glutardialdehyde, quaternary ammonium compounds.
- After using disinfectants, remove the disinfectant residue by wiping with a damp cloth.
- The surfaces must be dried immediately after disinfecting.

16.1.3 Removal of radioactive contaminants

- The agent must be specifically labelled as being an agent for removing radioactive contaminants.
- Ingredients of suitable agents for removing radioactive contaminants: anionic tensides, non-ionic tensides, polyhydrated ethanol.
- After removing the radioactive contaminants, remove the agent residue by wiping with a damp cloth.
- The surfaces must be dried directly after removing the radioactive contaminants.

16.2 Cleaning and disinfecting the haematocrit rotor

If despite all precautions a leakage of haematocrit capillaries during centrifugation should occur, rotor and lid must be cleaned according to the following instructions:

1. Carefully remove padded shells or sealing ring from the haematocrit rotor.
2. Soak haematocrit rotor, lid and padded shells or sealing ring in cold water, until the blood residues are completely dissolved.
3. Perform a suitable disinfection. The application of disinfection substances on the lid is to be avoided as the inscription can dissolve.
Ingredients of suitable disinfectants: glutaraldehyde, propanol, ethyl hexanol, anionic tensides, corrosion inhibitors.
Ingredients of suitable agents for removing radioactive contaminants: anionic tensides, non-ionic tensides, polyhydrated ethanol.
4. Rinse the components in cold water and dry them.
5. They can now be replaced into the centrifuge.

16.3 Procedure in case of broken haematocrit capillaries

Risk of injury!


Best work with two pairs of gloves worn on top of each other. Wear safety goggles and mask, as glass splinters could be blown off!

1. Carefully remove lid of haematocrit rotor.
2. Remove larger pieces of broken capillaries using pincers.
3. Remove rotor. Carefully and slowly remove the affected padded shells or the sealing ring with pincers and exchange them.
4. Clean and disinfect as described above.

16.4 Rotor removal / installation

Removal of haematocrit rotor:

1. Open the lid and disconnect the centrifuge from the power supply.
2. Remove the rotor lid by pressing the press lock.
3. Unscrew the fixing nut with the supplied spanner.
4. Pull the haematocrit rotor off the motor shaft.

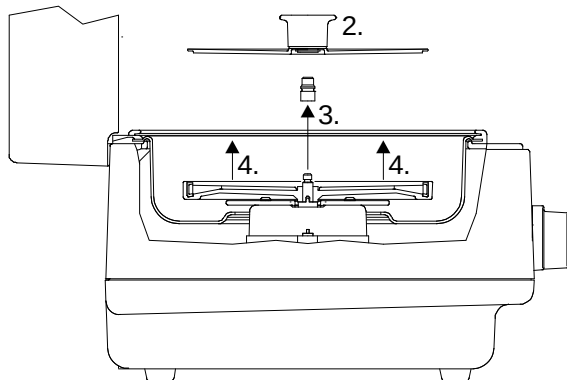Installation of the haematocrit rotor:

1. Insert the haematocrit rotor vertically onto the motor shaft. The motor shaft dog has to fit in the rotor slot.
2. Tighten the fixing nut with the supplied spanner.

The fixing nut must always be firmly tightened to the motor shaft.

17 Faults

If the fault cannot be eliminated with the help of the fault table, please inform Customer Service.

Please specify the type of centrifuge and the serial number. Both numbers can be found on the name plate of the centrifuge.

Perform a MAINS RESET:


- Switch off the mains switch (switch position "0").
- Wait at least 10 seconds and then switch on the mains switch again (switch position "I").

Fault	Display	Cause of fault	Remedy
No display	---	No voltage Mains input fuses defective.	- Check distribution voltage. - Check mains power input fuse, refer to Chapter "Change mains input fuse". - Mains switch ON.
System reset	- 2 -	Power failure during the centrifugation run. (The centrifugation run was not finished.)	- When stationary, open lid and press START key. - Repeat the centrifugation run if necessary.
Balance error	- 3 -	The rotor is unevenly loaded.	- Open lid when rotor is stationary. - Check the loading of the rotor, see chapter "Loading the rotor". - Repeat the centrifugation run.
Communication	- 4 -	Fault in control unit or power unit.	- Perform a MAINS RESET when the rotor has been stationary.
Overload	- 5 -	Motor or motor control defective.	- Perform a MAINS RESET when the rotor has been stationary.
Overvoltage Undervoltage	- 6 - - 8 -	Supply voltage outside tolerance (see Technical Data).	- Perform a MAINS RESET when the rotor has been stationary. - Check supply voltage.
Excess temperature	- 9 -	Excess temperature switch in motor has triggered.	- When rotor is stationary, open lid using emergency unlocking (see Emergency Unlocking chapter). - Allow motor to cool down.
Controller watchdog	- C -	Fault in control unit.	- Perform a MAINS RESET when the rotor has been stationary.
Lid error	- d -	Error in lid locking or lid closure.	- Perform a MAINS RESET when the rotor has been stationary.
Short circuit	- E -	Short circuit in control unit / power unit.	- Perform a MAINS RESET when the rotor has been stationary.
Configuration	- F -	Incorrect machine version	- Inform Customer Service.

18 Change mains input fuses

Switch off the mains switch and separate the device from the mains!

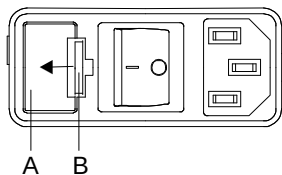The fuse holder (A) with the mains input fuses is located next to the mains switch.

- Remove the connecting cable from the machine plug socket.
- Press the snap-fit (B) against the fuse holder (A) and remove.
- Exchange defective mains input fuses.

Only use fuses with the rating defined for the type. See the following table.

- Reinsert the fuse holder until the snap-fit clicks shut.
- Reconnect the device to the mains supply.

Model	Type	Fuse	Order no.
HAEMATOKRIT 210	2104	T 3,15 AH/250V	E997
HAEMATOKRIT 210	2104-01	T 5 AH/250V	E914

19 Returning Devices

Before returning the device, a transport securing device has to be installed.

If the device or its accessories are returned to Andreas Hettich GmbH & Co. KG, in order to provide protection for people, the environment and materials, it has to be decontaminated and cleaned before being shipped.

We reserve the right to refuse contaminated devices or accessories.

Costs incurred for cleaning and disinfection are to be charged to the customer.

We ask for your understanding in this matter.

20 Disposal

Before disposal, the device must be decontaminated and cleaned to protect people, the environment and property.

When you are disposing of the device, the respective statutory rules must be observed.

Pursuant to guideline 2002/96/EC (WEEE), all devices supplied after August 13, 2005 may not be disposed as part of domestic waste. The device belongs to group 8 (medical devices) and is categorized in the business-to-business field.

The icon of the crossed-out trash can shows that the device may not be disposed as part of domestic waste.

The waste disposal guidelines of the individual EC countries might vary. If necessary, contact your supplier.

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All rights reserved. No part of this publication may be reproduced without the written prior permission of the copyright owner.